

# UV Flash

MIDBROOK MEDICAL



## 60 seconds

UV Flash can disinfect a wide range of everyday objects and items that are difficult to clean, including:

- Stethoscopes
- Oximeters
- Penlights
- PDAs
- Writing instruments
- Thermometers
- Cell phones
- Scissors
- Laptops
- Eyeglasses
- Blood pressure cuffs
- MARs scanners
- Oscopes
- Keys/ID cards
- Briefcases/satchels

## UVC light

At a certain wavelength, 254 nm, UVC breaks the molecular bonds inside microbial DNA, producing deadly thymine dimers, essentially, destroying the germ.

## Eliminate:

E.coli, MRSA, H1N1, SARS, pseudomonas aeruginosa, and many other pathogenic microorganisms. In studies going back more than 70 years, ultraviolet germicidal light (UVC) like that used in UV Flash® has been shown to kill more than 300 different germs, including bacteria, viruses, mold, fungi, and yeast. **There is no microorganism on earth that is resistant to ultraviolet germicidal light.**

UVC based infection prevention solution helps eliminate transient bacteria found on the everyday items and objects used and handled by healthcare workers. UV Flash is **designed to quickly and effectively disinfect** a diverse range of patient-care equipment as well as personal items of healthcare workers, all known to pose a risk of disease transmission and act as vector for hospital-acquired infections.

## Needed in:

- Hospital nursing stations
- ICUs
- Clinics
- Medical Offices, & More



**Disinfect and Sanitize with Light**

## UV Flash by Midbrook Medical

60 SECOND Infection Control Station

Phone: 1(517) 787-3481 • Email: Sales@midbrook.com

[www.midbrookmedical.com/uvflash](http://www.midbrookmedical.com/uvflash)


# UV Flash

MIDBROOK MEDICAL



## Needed in:

- Hospital nursing stations
- ICUs
- Clinics
- Medical Offices, & More

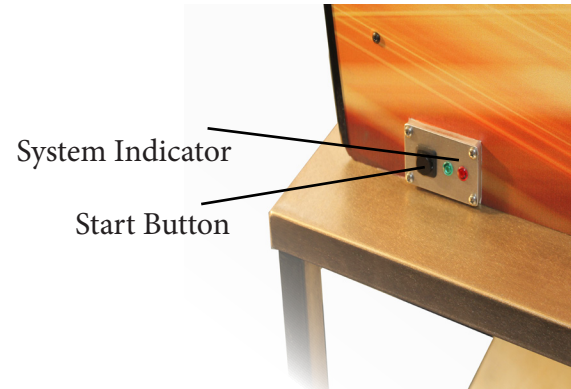

**“Kill E. coli, SARS, MRSA, in a FLASH.”**

## Review of UV Flash System Technical Specifications

|                         | System Only  | System with Cart |
|-------------------------|--|------------------|
| Dimensions/Weight       |  |                  |
| <b>HEIGHT</b>           | 36 in / 91.4 cm  | 70 in / 177.8 cm |
| <b>WIDTH</b>            | 22 in / 55.9 cm  | 28 in / 71.1 cm  |
| <b>DEPTH</b>            | 28 in / 71.1 cm  | 32 / 81.3 cm     |
| <b>WEIGHT</b>           | 60 lb / 27.3 kg  | 100 lb / 27.3 kg |
| Installation            |  |                  |
| <b>PLACEMENT</b>        | Cart-mounted on custom, stainless steel surgical case cart with antimicrobial copper alloy push handle (recommended) or countertop |                  |
| Electrical Requirements |  |                  |
| <b>CONNECTIONS</b>      | 120-volt, 20-amp circuit   |                  |
| Service Requirements    |  |                  |
| <b>CLEARANCE</b>        | UV Flash requires a minimum access of 2 inches (5cm) above the top   |                  |

## UV Flash by Midbrook Medical

60 SECOND Infection Control Station

Phone: 1(517) 787-3481 • Email: Sales@midbrook.com

[www.midbrookmedical.com/uvflash](http://www.midbrookmedical.com/uvflash)