


One Instrument - Three Configurations

Which Analyzer is Just Right for Your Laboratory?



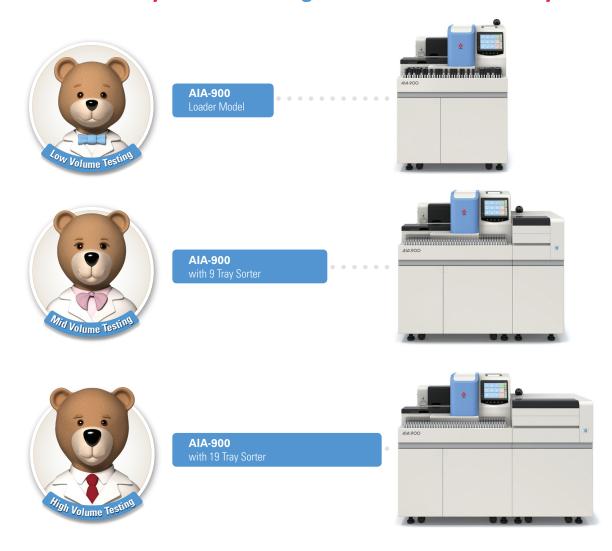

The latest addition to Tosoh's family of AIA immunoassay analyzers is available in three configurations, allowing you to choose a system that's *just right* for your laboratory workload.

- Instrument can be re-configured as testing volumes increase
- Throughput of 90 tests per hour
- Utilizes Tosoh's Unit Dose Test Cup reagent technology



TOSOH BIOSCIENCE

A System that Grows as Your Testing Volume Increases









AIA-900 Loader Model



AIA-900 with 9 Tray Sorter




AIA-900 with 19 Tray Sorter

Unit Dose Test Cup Technology

- Single, unitized cups use a dry reagent format that requires no pre-mixing, no pre-measuring, no on-board refrigeration and no waste.
- Dry reagent format ensures calibration stability of up to 90 days for most assays.
- · Reagents have a one year shelf life from date of manufacture.
- The AIA test cups and trays are labeled with code and lot number for automated scheduling and inventory management on AIA-900 sorter models.
- ST AIA-PACK reagents provide a faster turnaround time due to a 10 minute incubation protocol.
- Extensive test menu* provides a wide range of immunoassay tests including thyroid, cardiac
 markers, tumor markers, reproductive hormones, anemia markers, metabolic markers, diabetes,
 kidney markers and more.



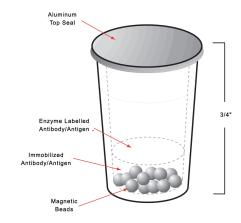

Test Cup Diagram



TOSOH BIOSCIENCE

Tosoh Bioscience, Inc. 6000 Shoreline Court, Suite 101 South San Francisco, CA 94080 Tel: (800) 248-6764 Fax: (650) 615-0415

www.tosohbioscience.us